

Anton Paar Kaomi for Nova

version 1.05
Copyright ©2020 Anton Paar GmbH



Report date:
File Name:

12/01/2025
BRF350 01-12-2025.qcuPhysIso

Operator:

labuser

Multipoint BET Summary/Results


Isotherm Branch Adsorption
Correlation coeff., r 0.999916

Slope 3194.12
C constant 46.2429

Intercept 70.5994
Surface area 1.067 m²/g

BET-Multipoint BET

—+— BET (All points) —○— BET function — Best Fit

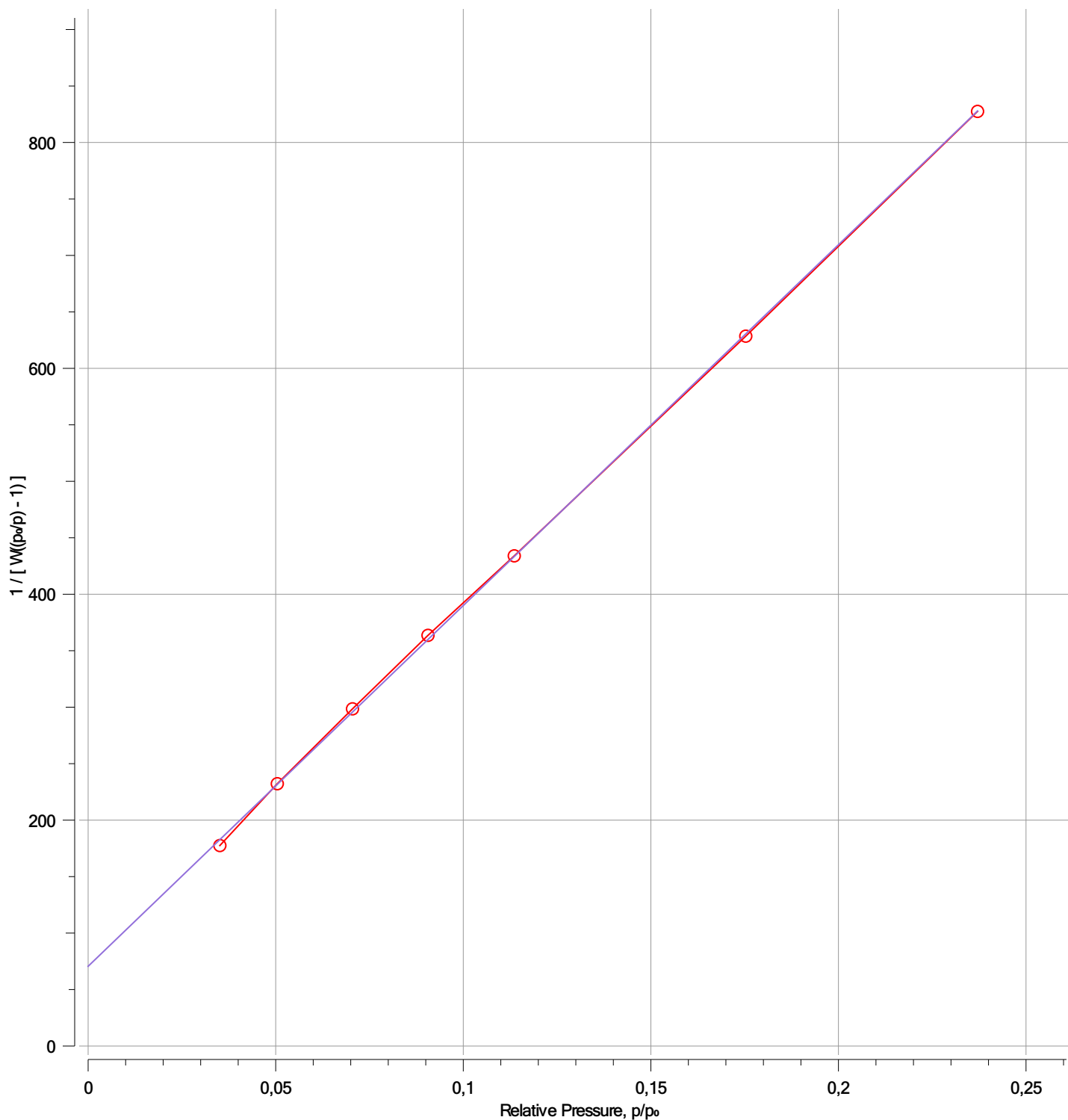

Anton Paar Kaomi for Nova

version 1.05
Copyright ©2020 Anton Paar GmbH



Report date: 12/01/2025 Operator: labuser
File Name: BRF350 01-12-2025.qcuPhysIso

BET-Multipoint BET		
Relative Pressure, p/p ₀	Volume Adsorbed @STP cm ³ /g	1 / [W((p ₀ /p) - 1)]
0.0351144	0.164035	177.5096
0.0504259	0.182951	232.2414
0.0704733	0.203191	298.5448
0.0905939	0.219209	363.6077
0.113591	0.236254	433.9894
0.175301	0.270598	628.5137
0.237084	0.300470	827.5110

Anton Paar Kaomi for Nova

version 1.05
Copyright ©2020 Anton Paar GmbH

Report date: 12/01/2025 Operator: labuser
File Name: BRF350 01-12-2025.qcuPhysIso

Analysis Data

Sample

ID	brf350 28-11-2025	Weight	1.6155 g
Name	brf350 28-11-2025	Description	brf350 28-11-2025

Analysis

Data ID	{fd157218-1897-433d-940c-f240b43803e8}		
Analysis Profile	N2 10 PT BET (biochar)		
Operator	labuser	Date	12/01/2025
Instrument	St 1 on NOVA 800 [s/n:1050059864]	Duration	170.52 min
Ambient Temp.	27.01 °C	Firmware	1.05
Warm Zone	3.0987 cm³	Void Volume Mode	He Measure
Thermal Delay	180 sec	Cell Type	9 mm with filler rod
		p ₀ Mode	From Ambient Pressure
		Cold Zone	22.1266 cm³

Adsorbate

Name	Nitrogen	Molecular Weight	28.0134 g/mol	Cross Sectional Area	16.2 Å²/molecule
Non-Ideality	6.58e-05 1/Torr	Bath Temperature	77.35 K		

Degas information

Type	Vacuum Degassing
Operator	labuser
Description	brf350 28-11-2025
Heating	Heat to 150.0 °C at 10.0 °C/min then hold for 600 min

Data Reduction Parameters

Data Reduction Parameters

Thermal Transpiration no

Adsorbate Model

Name	Nitrogen	Molecular Weight	28.0134 g/mol	Cross Sectional Area	16.2 Å²/molecule
Bath Temperature	77.35 K				